

Catalog #: U2269

Aliases

STING1; Stimulator Of Interferon Response CGAMP Interactor 1; STING; ERIS; MITA; Transmembrane Protein 173; TMEM173; NET23; MPYS; Endoplasmic Reticulum Interferon Stimulator; Stimulator Of Interferon Genes Protein; Endoplasmic Reticulum IFN Stimulator; FLJ38577; HSTING; HMITA; N-Terminal Methionine-Proline-Tyrosine-Serine Plasma Membrane Tetraspanner; Stimulator Of Interferon Response CGAMP Interactor-DeltaC; Stimulator Of Interferon Response CGAMP Interactor-DeltaN; Mitochondrial Mediator Of IRF3 Activation; Stimulator Of Interferon Protein; Stimulator Of Interferon Genes; Mediator Of IRF3 Activation; STING-Beta; Sting 1; SAVI

Background

Gene Name: STING1

NCBI Gene Entry: 340061

UniProt Entry: Q86WV6

Application Information

Molecular Weight: Predicted, 42 kDa; observed, 38 kDa

Clonality: Rabbit polyclonal antibody

Species Reactivity: Human, mouse, rat

Applications Tested: Western blotting (WB)

Immunogen

A synthesized peptide derived from human Phospho-STING (Ser366)

Isotype

Rabbit IgG

Storage Buffer

Supplied in PBS containing 50% glycerol, 0.5% BSA and 0.02% sodium azide.

Storage


Store at -20 °C for one year.

Recommended Dilutions

Western Blotting (WB): 1:500-1:2,500

Note: This product is for research use only.

Validation Data

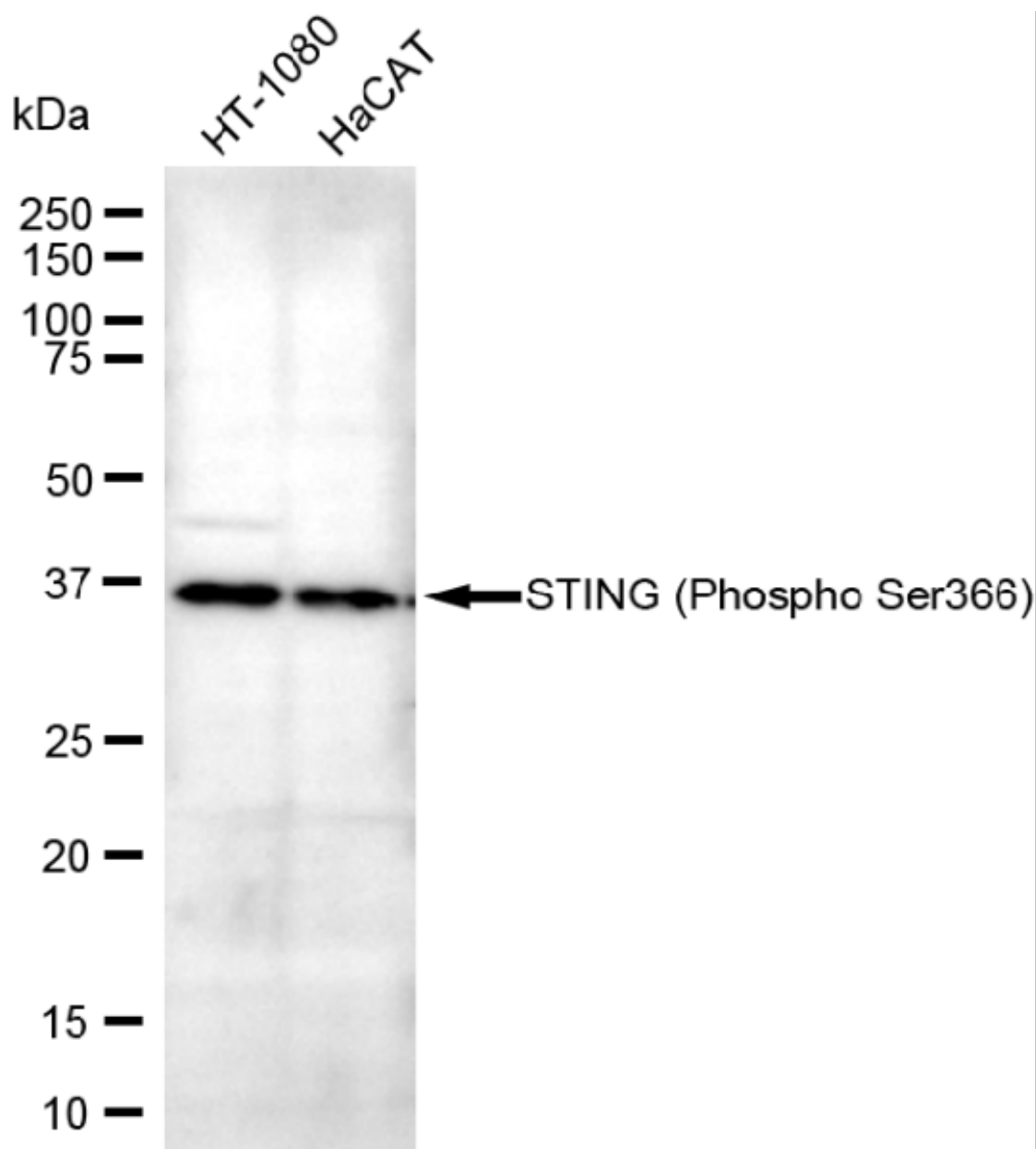

Copyright ©2025 Genuin Biotechnologies LLC

Western blotting analysis using anti-phospho-STING (S366) antibody (Cat#U2269). Total lysates (30 µg) were loaded and separated by SDS-PAGE. The blot was incubated with anti-phospho-STING (S366) antibody (Cat#U2269, 1:2,500) and HRP-conjugated goat anti-rabbit secondary antibody (Cat#201, 1:20,000) respectively. Image was developed using FeQ™ ECL Substrate Kit (Cat#226).

SUPPORT

SUPPORT@GENUINBIOTECH.COM
TEL: +1-540-855-7041

ORDERS

SALES@GENUINBIOTECH.COM
FAX: +1-540-855-7041

WWW.GENUINBIOTECH.COM